

Evelyn Metzger¹, Jason W Reeves¹, Emily Killingbeck¹, Yan Liang¹, Mona Khan², Laura Van Gerven³, Joseph Beechem¹, Peter Mommaerts²
¹ NanoString Technologies Inc., Seattle, WA, USA. ² Max Planck Research Unit for Neurogenetics, Frankfurt, Germany. ³ Department of Neurosciences & Department of Microbiology, Immunology and Transplantation, KU Leuven, Leuven, Belgium; Department of Otorhinolaryngology, University Hospitals Leuven, Leuven, Belgium


Abstract

Anosmia is a common symptom of COVID-19 and often persists well after the acute phase of the disease. The loss of smell is thought to be the result of the effects of SARS-CoV-2 on cell types that underly olfactory function. These cell types—which include sustentacular cells and olfactory sensory neurons—exist in the olfactory mucosa, which consists of an archipelago of islands surrounded by respiratory mucosa. Our earlier work using the GeoMx[®] Digital Spatial Profiler characterized whole-transcriptome effects of SARS-CoV-2 in the olfactory epithelium of a postmortem case at the spatial resolution of hundreds of cells. The transcriptional dysregulation of sustentacular cells we identified raises the question of which genes in sustentacular cells are altered as a result of viral load.

We have now addressed this spatial single-cell question using the CosMx[™] Spatial Molecular Imager (SMI). We built a spatially-informed atlas from healthy and disease states. Our panel consisted of 984 host targets and 9 probes for SARS-CoV-2. In total, we measured 63,589,058 spatial transcripts in 401,233 cells. We classified cells into known cell types using a combination of seed profiles from publicly available scRNA-seq data, semi-supervised clustering, and visual inspection of cell classifications and comparison with histological expectations. The semi-supervised clustering algorithm allowed us to discern more nuanced cell types of the epithelium (olfactory vs. respiratory horizontal basal cells) and to identify cell types that were not classified in the reference data (such as suprabasal cells). Since SARS-CoV-2 infection results in degradation of host mRNAs in the host cells, our approach was flexible enough to capture heavily infected cell types not adequately reflected in the reference (e.g., infected secretory cells and ciliated cells). Following classification, we focused on 725 sustentacular cells, grouped them into virus negative (535), low viral load (121), and high load (69), and found ~120 differentially expressed genes (DEGs). Differences between these groups contain “classic” DEGs such as TMPRSS2 and genes related to inflammatory or myeloid signaling.

We conclude by discussing the role of sustentacular dysregulation in anosmia. While our atlas was able to answer the specific biological question herein, we are addressing a myriad of other scientific inquiries related to viral infection and loss of smell. In total, our results underscore the importance of building spatially-enriched transcriptomic atlases of both healthy and diseased states.

Profiling key regions for SARS-CoV-2 Infection

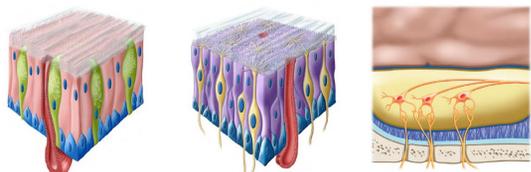

Respiratory mucosa (n = 1 control, 3 SARS-CoV-2+)
Olfactory mucosa (n = 2 control, 4 SARS-CoV-2+)
Olfactory bulb (n = 1 intact near olfactory mucosa)

Figure 1: Study design profiling various regions throughout the nasal mucosa
 Ten patients (11 tissue samples) were profiled using 984 host targets and 9 probes for SARS-CoV-2, including specific probes for the Delta and Omicron variants. These were sourced from both uninfected deceased patients (3) and patients with history of and evidence of SARS-CoV-2 infection (7) at the time of death. Below is the QC table showing the profiling depth and characteristics of each class of patient samples.

	Controls	Respiratory Mucosa	Olfactory Mucosa
Patients	3	3	4
Viral variant	-	2 omicron, 1 delta	2 omicron, 2 delta, 1 non VoC
Total Tissue Area (mm ²)	10.3	18.63	21.45
Number of Cells Analyzed	85,426	176,455	212,061
Total Transcripts Assigned to Cells	11.7M	28.7M	33.4M
Mean Transcripts per Cell	136.8	162.9	157.4
Maximum Transcripts per Cell	1,523	2,346	2,053

GeoMx identified preferential infection of sustentacular cells

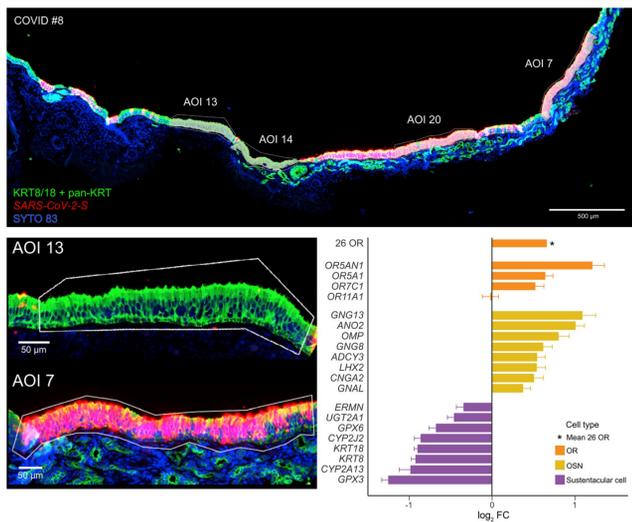


Figure 2: Transcriptional loss of sustentacular cell markers during SARS-CoV-2 Infection
 Profiling of highly infected SARS-CoV-2 patient with RNAscope & GeoMx Whole Transcriptome Atlas (WTA) in Khan et al [1] demonstrated preferential depletion of sustentacular cell markers from regions that were highly infected with the virus. Concurrent loss of expression of odorant receptor genes or markers of olfactory neurons was not observed.

Mapping interactions between SARS-CoV-2 infected cells and their surrounding environment

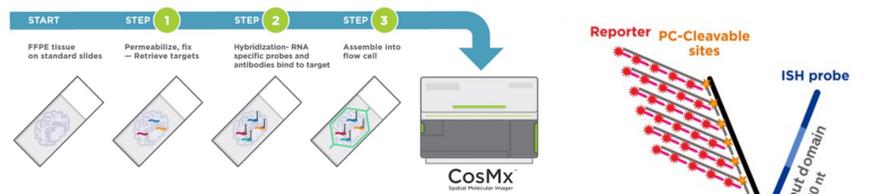


Figure 3: CosMx SMI profiling of FFPE samples
 Sequential flow of target-binding and fluorescently labeled reporter probes across FFPE tissues enable robust detection of RNA target fragments with <50 nucleotides of intact sequence required for target detection on standard histopathology slides.

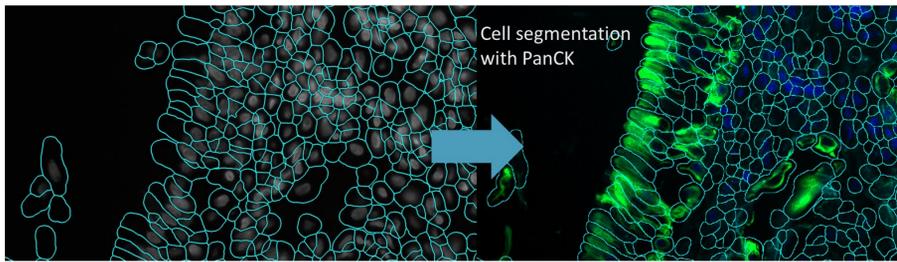


Figure 4 (Above): Segmentation of nasal epithelia for use with CosMx SMI workflow
 AI driven segmentation of epithelial cells performed with or without PanCK demonstrates improvements in segmentation of sustentacular cells.

Figure 6: Mapping respiratory & olfactory epithelial cell types in uninfected and infected deceased patients

Left) Common cell types making up the respiratory epithelium including ciliated cells, secretory cells, suprabasal cells, as well as cells populating the lamina propria of the nasal cavity.

Right) SARS-CoV-2-N (yellow points) expression localizes to infected cells making up the epithelia of the olfactory cleft preferentially. Infected cells primarily are sustentacular and ciliated cells

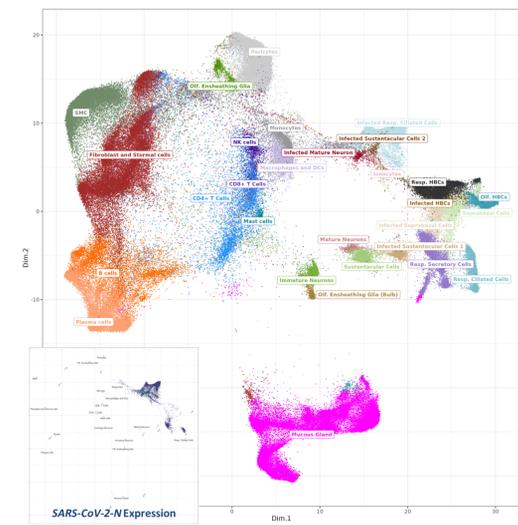
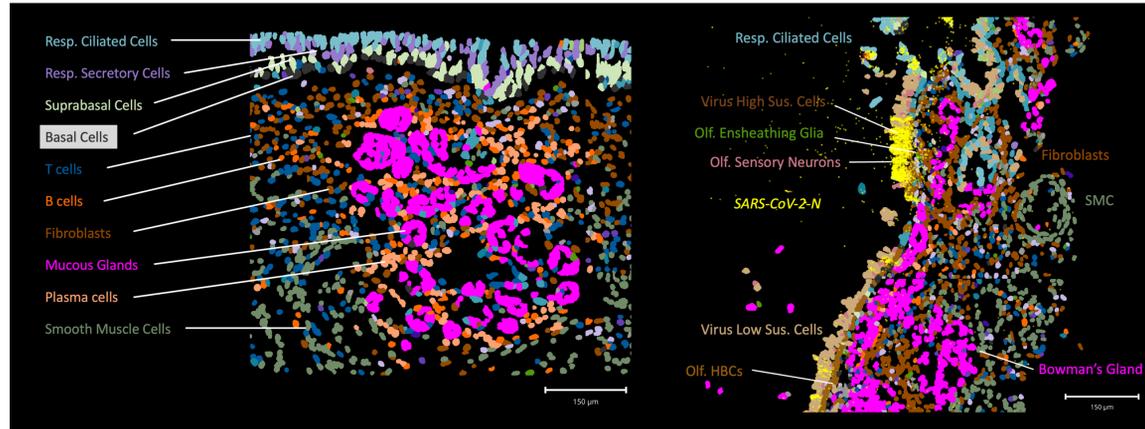


Figure 5: Atlas of SARS-CoV-2 infection in the respiratory & olfactory mucosa
 UMAP shows nasal & olfactory epithelia cell types identified in Durante et al [2] mapped to 401,233 cells profiled with CosMx SMI using InSituType algorithm [3]. Inset highlights infection of cells by SARS-CoV-2.

Anatomical barriers to SARS-CoV-2 transmission within the nasal cavity


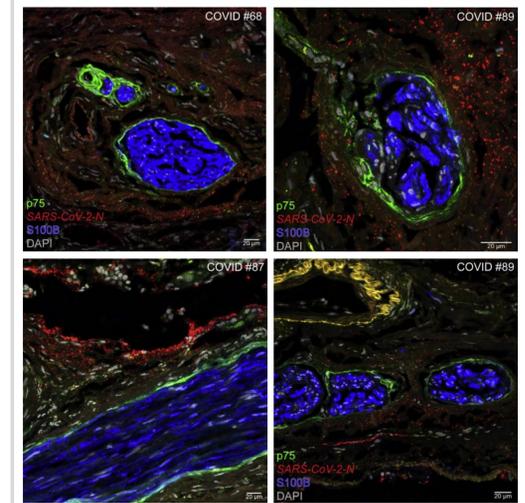


Figure 10: Anatomical barrier against SARS-CoV-2 neuroinvasion in the lamina propria of the olfactory mucosa

Further exploration with multiplexed IF & RNA Scope (Khan et al [4]) identified a novel set of fibroblasts (p75, green), associated with axon bundles (S100B, blue) in the lamina propria of the olfactory mucosa. These appear to prevent invasion of the virus (SARS-CoV-2-N red) by creating a physical barrier.

Figure 11: Identification of a novel type of barrier cells enwrapping olfactory axon fascicles of the olfactory mucosa

Cells surrounding the axon fascicles can be segregated out using CosMx SMI validating previous results suggesting their importance as well as robust characterization of their function *in situ*. Left) All cells, Right) Fibroblasts only + PanCK Immunofluorescence.

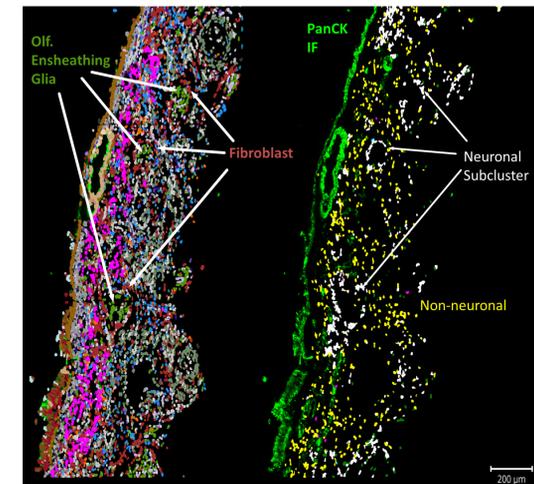

Conclusions

- CosMx SMI deeply characterized FFPE nasal epithelia from postmortem tissue samples using both the standard universal CosMx RNA panel and spike-ins for odorant receptor genes and viral genes
- Spatial profiling of nasal epithelium identifies local modifiers of SARS-CoV-2 infection and invasion
- Spatial cell typing identifies additional cell types not previously characterized using single-cell dissociated sequencing and consequences of viral infection
- 120 differentially expressed genes were identified between infected and uninfected sustentacular cells
- A barrier to neuro-invasion was identified specifically by spatial profiling

Scan here to download or learn more

The CosMx[™] SMI and decoder probes are not offered and/or delivered to the following UPC member states* for use in these countries for the detection of RNA in a method used for the detection of a plurality of analytes in a cell or tissue sample without the consent of the President and Fellows of Harvard College (Harvard Corporation) as owner of the Unitary Patent EP 4 108 782 B1. The use for the detection of RNA is prohibited without the consent of the of the President and Fellows of Harvard College (Harvard Corporation).
 *Austria, Belgium, Bulgaria, Denmark, Estonia, Finland, France, Germany, Italy, Latvia, Lithuania, Luxembourg, Malta, Netherlands, Portugal, Slovenia, Sweden

The CosMx[™] SMI and decoder probes are not offered and/or delivered to the Federal Republic of Germany for use in the Federal Republic of Germany for the detection of cellular RNA, messenger RNA, microRNA, ribosomal RNA and any combinations thereof in a method used in fluorescence in situ hybridization for detecting a plurality of analytes in a sample without the consent of the President and Fellows of Harvard College (Harvard Corporation) as owner of the German part of EP 2 794 928 B1. The use for the detection of cellular RNA, messenger RNA, microRNA, ribosomal RNA and any combinations thereof is prohibited without the consent of the of the President and Fellows of Harvard College (Harvard Corporation).