

Vantage 3D RNA:Protein Solid Tumor Assay for FFPE

The nCounter® Vantage 3D RNA:Protein Solid Tumor Assay for FFPE simplifies RNA and protein expression analysis with curated content covering up to 770 RNA targets and 26 total and phospho-protein targets. This highly-multiplexed assay is capable of simultaneously characterizing RNA and protein expression from as little as two FFPE slides.

The core nCounter technology uses unique molecular barcodes to detect nucleic acids of increasing variety. Specifically, antibodies of interest are barcoded with unique synthetic DNA oligonucleotides with a photocleavable linker. Each DNA oligonucleotide is then UV cleaved and recognized by a unique Reporter probe that contains a fluorescent barcode. The fluorescent probes are then imaged and counted by the nCounter Analysis System to provide a direct, digital readout of protein expression. The result is an integrated RNA:Protein workflow.

Learn more about [3D Biology™ Technology](#).

Product Workflow

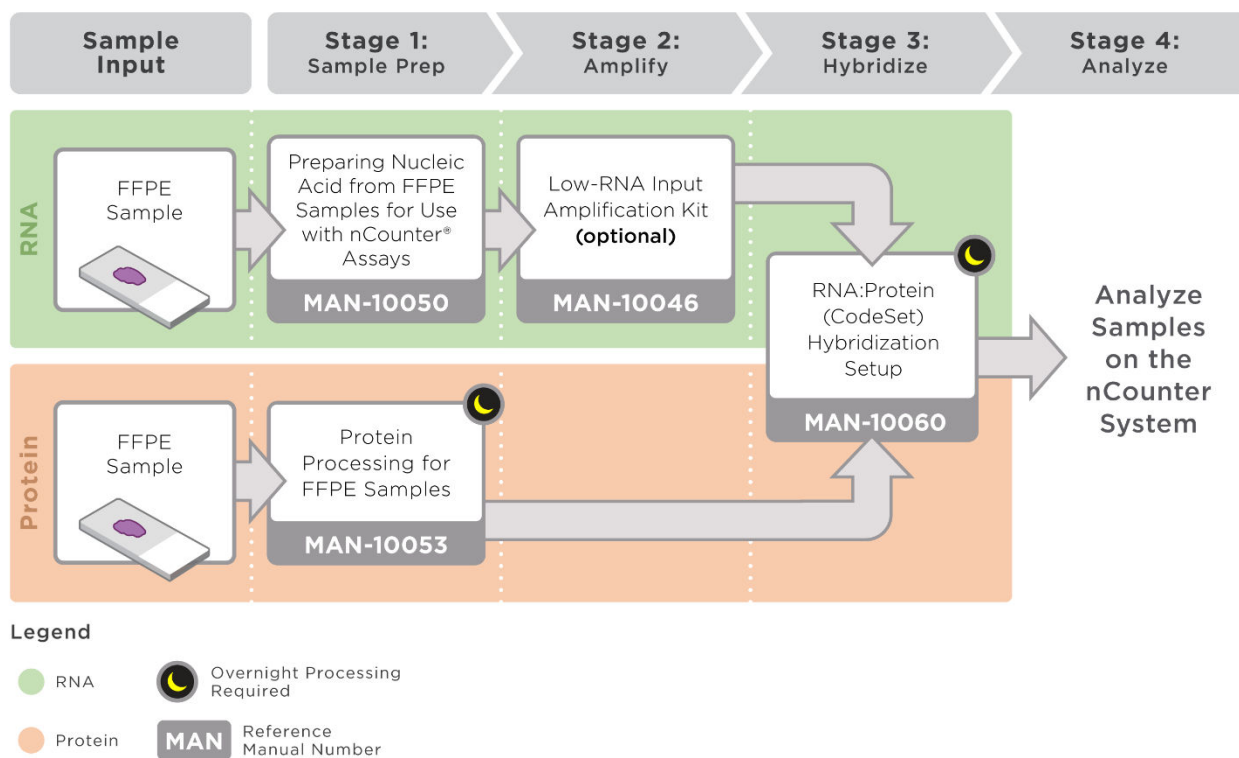


Figure 1. Workflow for Vantage 3D RNA:Protein Solid Tumor Assay for FFPE

FOR RESEARCH USE ONLY. Not for use in diagnostic procedures

© 2016–2017 NanoString Technologies, Inc. All rights reserved.

NanoString, NanoString Technologies, the NanoString logo, nCounter, Vantage 3D, and 3D Biology are trademarks or registered trademarks of NanoString Technologies, Inc., in the United States and/or other countries.

Materials and Supporting Documents

Table 1. Materials provided in the Vantage 3D RNA:Protein Solid Tumor Assay for FFPE Kit

Kit	Reagents	Storage
Vantage 3D RNA:Protein Solid Tumor Assay for FFPE Catalog #: VRPC-SPKP-HSTF-12	RNA	
	Reporter CodeSet	-80°C
	Capture ProbeSet	-80°C
	Protein	
	Protein TagSet (R)	-80°C
	Antibody Mix	-80°C
	Buffer W	4°C
	Buffer T	4°C

NOTE: Please reference the manuals listed in [Figure 1](#) and [Table 2](#) for additional required reagents not supplied by NanoString.

Table 2. Supporting Documents

Step	Manual	Protocol
Nucleic Acid Extraction	MAN-10050	Preparing Nucleic Acid from FFPE Samples for Use with nCounter Assays
RNA Amplification (optional)	MAN-10046	Low-RNA Amplification Kit
Protein Preparation	MAN-10053	Protein Processing for FFPE Samples
Hybridization	MAN-10060	RNA:Protein Hybridization Setup (CodeSet)

Intellectual Property Rights

This nCounter Vantage 3D RNA:Protein Solid Tumor Panel for FFPE Overview and its contents are the property of NanoString Technologies, Inc. ("NanoString"), and are intended for the use of NanoString customers solely in connection with their operation of the nCounter Analysis System. The nCounter Analysis System (including both its software and hardware components) and this User Manual and any other documentation provided to you by NanoString in connection therewith are subject to patents, copyright, trade secret rights, and other intellectual property rights owned by or licensed to NanoString. No part of the software or hardware may be reproduced, transmitted, transcribed, stored in a retrieval system, or translated into other languages without the prior written consent of NanoString. For a list of applicable patents, see www.nanostring.com/company/patents.

Limited License

Subject to the terms and conditions of sale of the nCounter Analysis System, NanoString grants you a limited, non-exclusive, non-transferable, non-sublicensable, research use only license to use this proprietary nSolver™ software with the nCounter Analysis System only in accordance with this manual, the manual for the nCounter Analysis System, and other written instructions provided by NanoString. Except as expressly set forth in the terms and conditions, no right or license, whether express, implied, or statutory, is granted by NanoString under any intellectual property right owned by or licensed to NanoString by virtue of the supply of this software or the proprietary nCounter Analysis System. Without limiting the foregoing, no right or license, whether express, implied, or statutory, is granted by NanoString to use the nSolver Analysis Software or nCounter Analysis System with any third-party product not supplied or licensed to you by NanoString, or recommended for use by NanoString in a manual or other written instruction provided by NanoString.